

Sterilization Tubes & Pouches

Chex-All® Sterilization Pouches and Tubes are available in a variety of sizes and used to contain small or individual instruments and implants during steam or ethylene oxide sterilization.



All varieties feature:

- Triple sealed paper/plastic film: Provides excellent sterility protection
 - Premium quality 63.5µ plastic and 73-74g medical grade paper
 - Available with heat sealing or self sealing closure
 - Wide adhesive line guarantees easy and reliable sealing
 - Chevron with corner seals assures dust-free, aseptic presentation
 - For Steam and Ethylene Oxide Sterilization
- Steam:** Yellow to grey/ black color change
Ethylene Oxide: Red to green color change

Chex-All® Heat Sealing Tubes & Pouches

Premium quality sterile packaging products with independent indicators imprinted on both inside and outside of every pouch. External indicators serve as process indicators. Internal indicators monitor the effective penetration of steam or EO gas on the inside of the package.

Each roll of Chex-All® Tubes measures 100ft (30.5m) per roll.

Chex-All Pouches

Reorder no.	Qty. per case	Inches	Centimetres
02500800	1000	3x8	7.6x20.3
02501000	1000	5x10	12.7x25.4
02501200	1000	7x12	17.8x30.5
02501400	500	8x16	20.3x40.6
02501800	500	12x18	30.5x45.7
02502400	500	18x24	45.7x60.9

Chex-All Tubes

Reorder no.	Qty. per case	Inches	Centimetres
02600200	10	2	5.1
02600300	10	3	7.6
02600400	10	4	10.2
02600600	10	6	15.2
02600900	6	9	22.9
02601200	6	12	30.5



Sterilization Tubes & Pouches


Chex-All® II Self-Sealing Pouches

Chex-All® II Pouches contain internal and external chemical indicators, that experience a yellow to grey/black color change in steam and red to green color change with the presence of ethylene oxide. External indicators demonstrate that the pouch has been exposed to steam or ethylene oxide gas. Internal indicators ensure the sterilant has penetrated the pouch and reached the contents.

Chex-All® II Self Sealing Pouches save time and are easy to use!



Reorder no.	Quantity per case*	Inches	Centimetres
02400800	500	3x8	7.6x20.3
02401000	500	5 x10	12.7x25.4
02401200	500	7x12	17.8x30.5
02401400	250	8x16	20.3x40.6
02401500	500	5x15	30.5x45.7
02401800	250	12x18	45.7x60.9

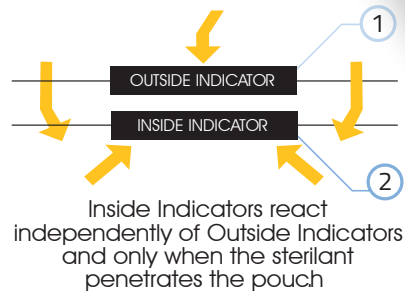

Chex-All® III Self-Sealing Pouches

Process indicators are printed on the pouch corners and covered by the extension of the transparent film beyond the sealed area. This guarantees protection of the indicators from direct contact with metal trays or instruments and enables clear visualization of results.

Pre-printed process indicators experience a easy-to-interpret color change; yellow to grey/black in steam and red to green in EO.

Chex-All® III pouches are packed in small boxes for easy distribution between departments



Reorder no.	Quantity per case	Inches	Centimetres
02430900	1000	3.5 x 9.5	8.9x24.1
02431000	1000	5.25 x 10	13.3x25.4
02431300	1000	7.5 x 13	19x33.02
02431400	500	9 x 16	20.3x40.6
02431800	500	13 x 18	33x45.72

