

Steri-Wrap® Unite 60

Sterilization Wrap

Product Description

Propper Steri-Wrap Unite 60 sterilization wrap is composed of two sheets of SMS bonded together for convenient simultaneous wrapping of one or a collection of medical devices that will be sterilized following standard healthcare practices. The wraps are made of polypropylene with blue or green pigment and an anti-static treatment. Steri-Wrap Unite is used to contain instruments during steam, ethylene oxide and hydrogen peroxide sterilization and maintain sterility after processing.

Steri-Wrap Unite 60 is available in sheets

31502400 Steri-Wrap® Unite 60 - 24" x 24" Wrap - 1 Case of 120
31503600 Steri-Wrap® Unite 60 - 36" x 36" Wrap - 1 Case of 72
31504800 Steri-Wrap® Unite 60 - 48" x 48" Wrap - 1 Case of 48
31505400 Steri-Wrap® Unite 60 - 54" x 54" Wrap - 1 Case of 48


Benefits

Steri-Wrap Unite 60 sterilization wraps allow effective sterilization penetration and maintain sterility of reusable devices during handling and storage. The two-color sheeted application allows for easier visualization of damage, tears or punctures and can be used to differentiate between loaner and facility owned trays.

Steri-Wrap Unite 60 conforms to ISO 11607-1:2019 and EN 868-2:2017

Technical Data

Property	Steri-Wrap® Unite 60	Unit
Dry Mullen Burst Strength	3059	g/cm ²
Dry Grab Tensile Strength CD	12000	grams
Trapezoid Tear MD	4000	grams
Hydrostatic Resistance (60cm/min)	70	cm
Basis Weight	60	gsm

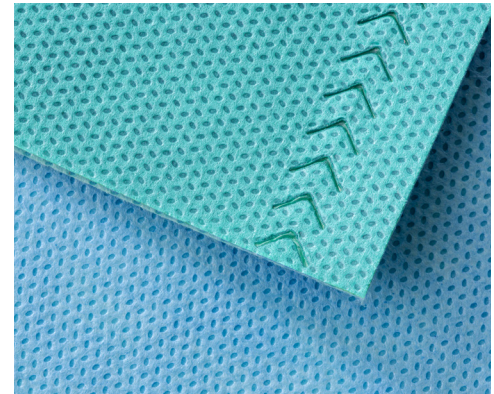

Shelf Life

60 months

Shelf Life After Sterilization

365 days

Storage Conditions

Cool dry place out of direct sunlight.



36-04 Skillman Avenue
Long Island City, NY 11101
Tel. (718) 392-6650 • Fax. (718) 482-8909
www.proppermfg.com