

Duo-Flash® Sterilization Indicator

Duo-Flash® is a 4 inch (10 cm) dual cycle chemical indicator for immediate use steam sterilization monitoring.

For reliability you can see and results you can count on.



Duo-Flash®

For use in Steam Sterilizers

The Duo-Flash® immediate use steam sterilization (IUSS) indicator is specifically engineered for two uses: nonporous loads and porous loads.


For nonporous loads: After monitoring a 3-minute, 272°F(133°C) IUSS cycle used for nonporous load, the arrow pointing to the 3 turns from green to dark.

For porous loads: After a 10-minute, 272°F(133°C) IUSS cycle used for a porous load, the arrow pointing to the 10 turns from green to dark.

Duo-Flash® has an extra measure of security

Plastic lamination on both sides maintains indicator integrity in demanding unwrapped cycles. Its durable construction, convenient design, and permanent color change were all developed with the requirements of IUSS in mind.

If required, the indicator may be retained for record keeping for as long as necessary.

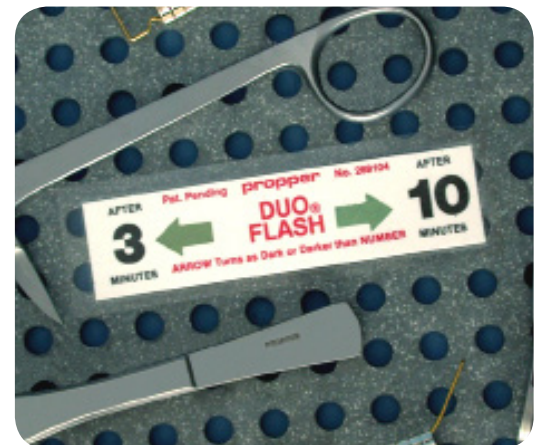

| Reorder No. I | Qty |
|---------------|-----|
| 26910400 | 100 |



Propper Manufacturing Company, Inc | marketing@proppermfg.com | www.proppermfg.com

36-04 Skillman Ave. Long Island City, NY, USA 11101 | Tel : +1-718-392-6650 | Fax : +1-718-482-8909