


# Vapor Line® Steam Sterilization Integrators

The Vapor Line® Type 5 Steam Integrator is an immediate readout chemical indicator for monitoring the essential conditions of steam sterilization in gravity & pre-vacuum sterilizers.

**Vapor Line®** is designed to react to all critical variables of steam sterilization and is equivalent in performance to biological indicators. Vapor Line® features an overkill security margin to provide added protection in infection prevention.

- During sterilization the dark bar irreversibly migrates along the Pass/Fail channel
- Any migration further than the Fail Area indicates sterilization conditions were achieved
- For the maximum protection, Vapor Line® is shipped in a foil pouch
- Vapor Line® may be used in a process challenge device as part of release criteria for nonimplant loads and for loads containing implants in an emergency

(AAMI ST-79:2017, ISO 11140-1, Type 5)

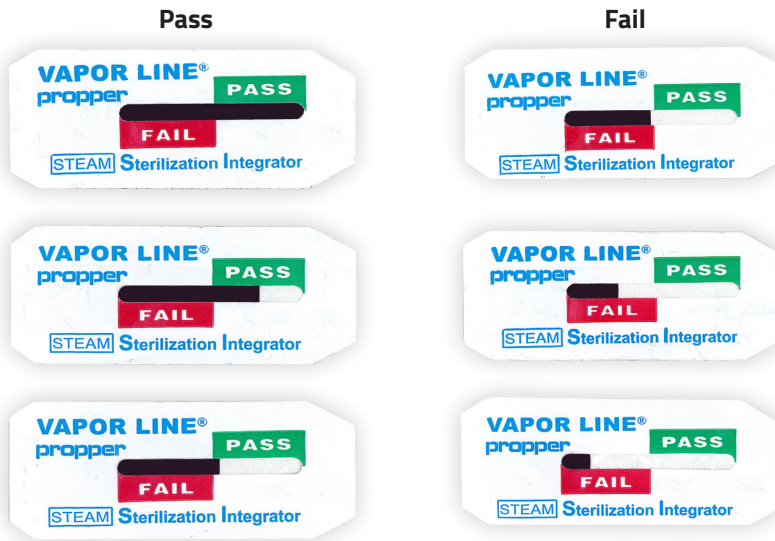

Reorder No	Description	Quantity
26900510	Without Extender 5.1cm x 2.2cm 2" x 7/8"	100 integrators/pouch
26900925	Without Extender 5.1cm x 2.2cm 2" x 7/8"	250 integrators/pouch
26900825	With Extender 22.9cm x 2.2cm 9" x 7/8"	250 integrators/pouch

## Vapor Line® Extender

The Extender elongates the Vapor Line® integrator for easier visualization in the tray and easier pack retrieval after processing

Vapor Line® Extenders have a fill-in area for recording batch-specific information, ideal for record keeping and a chemical indicator which exists a white to black color change when sterility conditions are met.



## VAPOR LINE® Extender

LOT NO.  
DATE



Steam Indicator

